BRUKSANVISNING / INSTRUCTION MANUAL

NAKKEMASSASJE / NECK MASSAGER / 80000913

INNLEDNING

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INTRODUCTION

Thank you for choosing a product from Epiq. You have acquired a quality product which meets all applicable performance and safety standards in the EU. In order to ensure proper treatment and a long product life, we recommend that you comply with the instructions below.

Please read the manual and especially the safety instructions carefully before using the product. Keep this manual and provide it for referral along with the product when others are using it.

GENERAL SAFETY INSTRUCTIONS

- Improper use of the product can cause personal injury and damage.
- Only use the product for what it is intended. The manufacturer is not liable for damages arising from improper use or handling.
- Always place the product on a dry, flat and stable surface well away from table edge or similar.
- The product is intended for household use only. The product is not for professional use.
- The product may only be connected to 220-240V~, 50/60Hz.
- The product engine/electronics, cord or plug cannot be submerged in water, and no water may get in contact with the engine/electronics.
- Never touch the product, cord or plug with wet or damp hands.
- Do not expose the product to direct sunlight, high temperatures, humidity, dust or corrosive substances.
- Never leave the product when it is turned on.
- When the product is in use, it should be kept under constant supervision. Children should always be supervised when using the product to ensure that they do not play with the product. The product is not a toy.
- The product may have sharp parts and must, therefore, be kept

out of reach of children.

- Only use the accessories that come with the product or that are recommended by the manufacturer.
- The product must not be used by people with reduced sensitivity, physical or mental disability, or people who are not able to operate the product unless they are supervised or instructed by a person responsible for their safety.

CORD AND PLUG

- Do not let the cord hang over the edge.
- Keep the cord and the product away from heat sources, hot objects and open flames.
- Unplug the power cord before cleaning and if the product is not in use. Do not pull the cord when removing the plug from the socket, rather grasp the plug itself.
- Make sure that the cord is fully extended.
- Check regularly that the cord and the plug are not damaged and do not use the product if this is the case, or if it has been dropped on the floor, dropped in water or damaged in another way.
- If the product or plug is damaged, the product must be inspected and, if necessary, repaired by an authorised repairer. Never try to repair the product yourself. Please, contact the store to see if there is a repair that falls under warranty. Tampering of the product will void the warranty.

RULES OF COMPLAINTS THE WARRANTY DOES NOT APPLY:

- If the above is not observed.
- If there has been unauthorised intervention of the product.
- If the product has been mishandled, subjected to rough treatment or some other form of damage.

PRODUCT DISPOSAL


Used product should not be disposed as household waste. According to legislation, the product should be disposed of

properly so that the materials can be recycled in an environmentally friendly manner. Further information is available from the responsible local authority or local disposal company. The product is equipped with the following logo:


PRODUCT ACCESSORIES

1. The main belt of massager
2. The control panel
3. The power cord
4. The power plug
5. The hand pocket

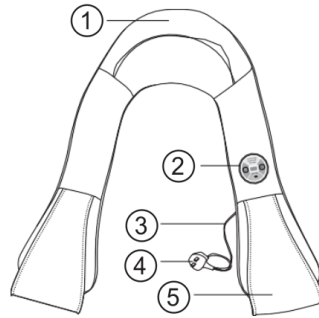THE CONTROL PANEL

1. Countdown/gear display
2. On/Off
3. Modes
4. The intensity decrease button
5. The intensity increase button

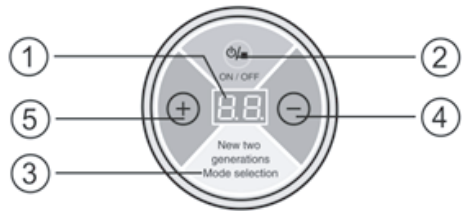PRODUCT FEATURES

- Two-point beating-massage can relieve shoulder and neck pain.
- Function: 99 kinds of people-like beating-massage modes, 20 different gears of intensity adjustment, automatic timing 10 minutes.

PRODUCT PERFORMANCE AND EFFECTS

- Relieve cervical pressure, eliminate and reduce cervical pain.
- Promote blood circulation, slow muscle pain and stiffness.
- Ease tension, reduce fatigue and burnout.
- Promote cell movement and improve the blood circulation.

THE OPERATION METHOD

1. The preparation before use

Insert the power plug into the power outlet, the digital display shows OFF, and the massager is in standby mode at this time, and then place the massager at the correct position. [Please refer to the "The Using Method"].

2. Automatic mode

Press [ON/OFF] to turn on the massager, then massage ribbon is started and enters into the automatic mode to beat.

3. Manual mode

Press the button to enter the manual mode, each time the button is pressed, the beating-massage modes changed accordingly.

After all the modes of the beating-massage is finished, it returns to the automatic state and this cyclic process is repeated.

4. The adjustment of the intensity

Press the button to increase the intensity of the beating-massage, while press the button to reduce the intensity, and at the same time, the corresponding intensity levels are displayed.

5. The termination of massage

The default massage time is ten minutes, when turning on the massager, it enters into the countdown state, and will be automatically shut down after ten minutes.

THE USING METHOD

Neck/Shoulder

Back

Back

Waist

Thigh

TECHNICAL DATA

Article No. 80000913
220-240V~, 50/60Hz
60 Watt

PRODUCER

Epiq ApS
Krøyer Kielbergs Vej 3, 3. th
8660 Skanderborg, Denmark
www.epiq.dk

We cannot be held responsible for any printing errors.